

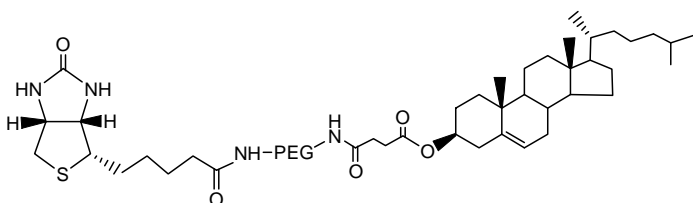
TECHNICAL DATA SHEET**Cholesterol PEG Biotin, MW 1k, 2k, 3400, 5k, 10k, 20k**

Catalog Numbers: PG2-BNCS-1k, 2k, 3k, 5k, 10k, 20k.

Synonym: Cholesteryl PEG Biotin, Biotin PEG Cholesterol

Description:

Cholesteryl PEG Biotin is one of Nanocs' biotinylated Cholesterol PEG derivatives with strong binding affinity to avidin or streptavidin. Cholesterol is an essential lipid component to build and maintain basic animal cell function. It is one of animal cell-membrane structure components and establishes proper membrane permeability and fluidity in animal cells. In addition to its structural function, cholesterol is also the important precursor for the biosynthesis of steroid hormones, bile acids and vitamin D. Cholesterol is a hydrophobic molecule with very low water solubility. Biotin functional Cholesterol PEG can be used for molecule detection, targeting and delivery through its specific interaction with avidin or streptavidin. **Cholesterol PEGs** offered by Nanocs have good water solubility, flexible linker structure and diverse functionalities. Our Cholesterol PEG derivatives have been used successfully in cell membrane imaging, liposome formulation and drug delivery. As a leading PEG cholesterol provider, Nanocs now offers a variety of multifunctional PEG cholesterol molecules for different applications.

Product Structure:**Product Specifications:**

- Composition: **Cholesterol PEG Biotin**
- Appearance: White/off-white solid, semi-solid depends on molecular weight.
- Solubility: 10 mg/mL, clear in water, chloroform, DMSO.
- Functional group: Cholesterol and Biotin

Handling and Use:

Cholesteryl PEG Biotin is sensitive to moisture and temperature. For best use, material should always be kept in low temperature in dry condition. Prepare fresh solution right before use. Avoid frequent thaw and freezing. For more information about using this product, visit www.nanocs.net.

Storage Conditions:

Biotin PEG Cholesterol should be stored at -20 °C. Desiccate. Materials may be handled under inert gas for best stability. Re-test material after 6 months.

This product is for research use only and is not intended for use in humans or for diagnostic use.

Related Products:

Cholesterol PEG NHS	Cholesterol PEG FITC
Cholesterol PEG DBCO	Cholesterol PEG Azide
Cholesterol PEG COOH	Cholesterol PEG Amine
Cholesterol PEG Maleimide	

To Order:Order online: www.nanocs.netOrder by Email: sales@nanocs.com

Order by phone: 1(800) 388-4221; 1(888) 908-8803

For more information, visit www.nanocs.net